

Air Scaler

TA-100

4H / 2H

The standard model air scaler that offers powerful and effective tartar removal.

Smart and powerful




Simple setup just to connect to your turbine hose makes TA-100 ready for powerful and effective scaling. The subsonic vibration provides high working efficiency and patient comfort. The compact and sleek design fits in hand comfortably and reduces hand fatigue. The slimmer scaling tips enter hard-to-access areas like subgingival easily.

Lightweight stainless steel body



Besides the main body, the joint of TA-100 is made of stainless steel. This is to provide TA-100 with higher resistance against strain that is applied to the joint when connecting and disconnecting. Robustness and lightweight are combined in TA-100. Less than 60 grams.

Powerful scaling

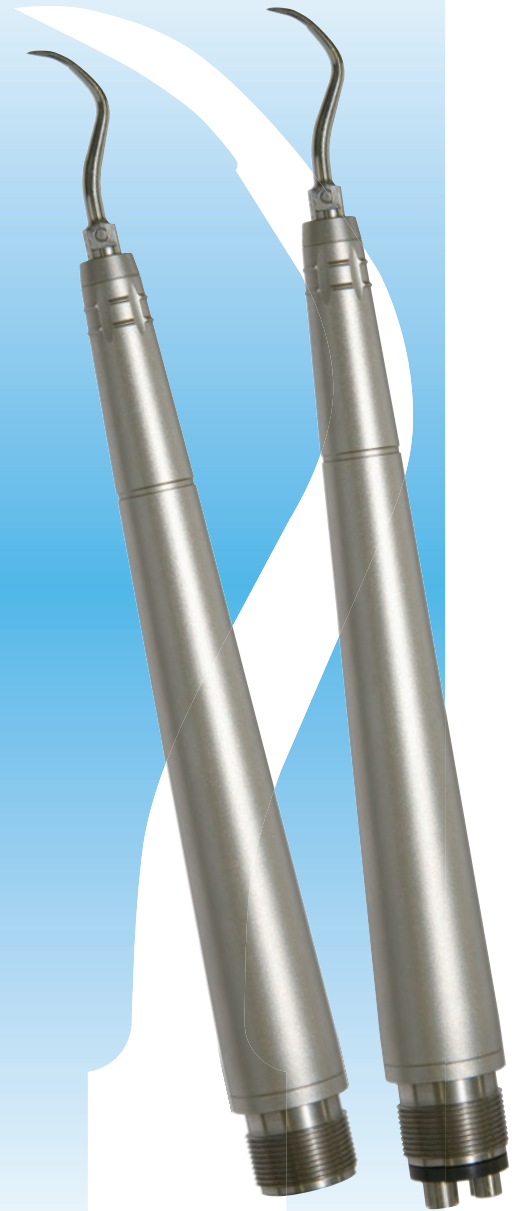

Specifications

Model	TA-100-4H	TA-100-2H
Backend configuration	MIDWEST 4 Hole	BORDEN 2 or 3 Hole
Frequency	about 6,000 Hz	
Coupling pressure	0.3 MPa (45 psi, 3.0 bar, 3.0 kgf / cm ²)	
Weight (without scaling tip)	58 g (2.1 oz)	59 g (2.1 oz)
Length (without scaling tip)	129 mm (5.1 inch)	123 mm (4.8 inch)

Specifications are subject to change without prior notice for improvements.

Manufacturer : MICRON CORPORATION

2-17-7, Ikegami, Ota-ku, Tokyo, 146-0082, Japan
 Tel : +81 3 3755 0396, Fax : +81 3 5747 5396
<http://www.micdent.com/english/> E-mail : office@micdent.com

European Authorized Representative : EMERGO EUROPE

Prinsessegracht 20 2514 AP The Hague The Netherlands
 Tel : +31 70 345 8570, Fax : +31 70 346 7299